


INTEGRATED SOLUTIONS

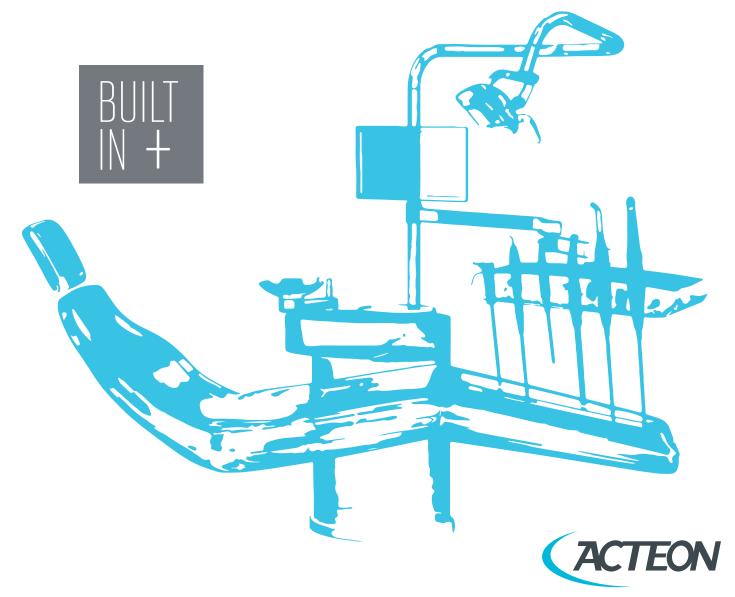

WHY INTEGRATE?

CLINICAL EXPERTISE

As a major partner for equipping dental chairs, ACTEON® places its technical and clinical expertise at the service of practitioners, from conception to delivery to be always closer to your needs.


FULL RANGE

ACTEON® provides a high-quality and affordable full range of innovative instruments.

LATEST TECHNOLOGY

Ultrasonics, polishers, curing lights, intraoral cameras and disposable air/water syringes.



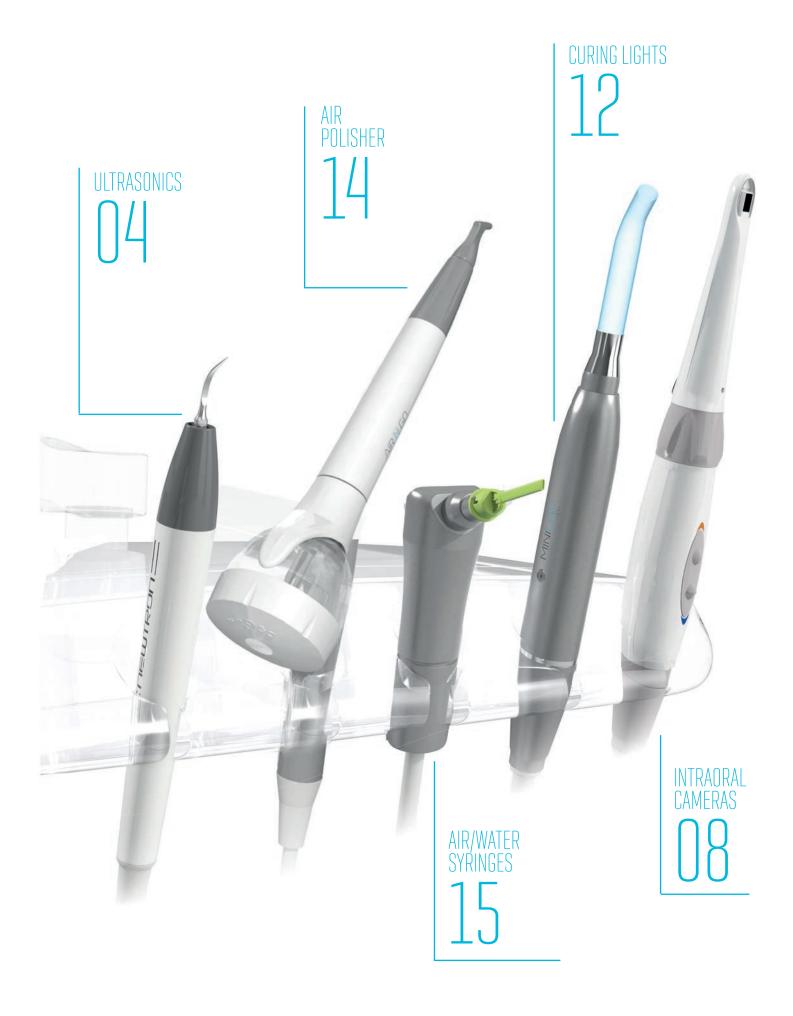

ULTRASONICS

NEWTRON MAKE THE SWITCH





COMFORT FOR YOU & YOUR PATIENTS

Very little pressure, with gentle vibrations

- less hand fatigue
- ideal for a wide variety of needs: from sensitive or disabled patients to pediatrics



FASTER TREATMENT TIME

- Treatments 20% to 50% shorter⁽¹⁾
- More time for hygiene education



DON'T GO BACK TO THE STACKS

Magnetostrictive inserts have fragile and tempermental stacks

- If bent or damaged, their power and performance decreases
- Not an issue with NEWTRON® tips or handpiece



LESS HEAT, LESS WATER

- Magnetostrictive stacks create more heat and need more water to cool handpiece, tip and tooth
- Decreased visibility for clinician
- · Can cause gagging for patients



(1) Hayes MJ, Cockrell D, Smith DR. A systematic review of musculoskeletal disorders among dental professionals. IFDH. 2009 Aug;(3):159–165




MORE VERSATILITY


NEWTRON® tips are gentle on teeth and have the most applications of any ultrasonic scaler

Over 80 different tips, with exclusive and patented designs, alloys and coatings that won't damage the different surfaces treated.



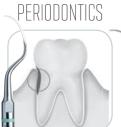



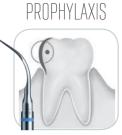





With over 40 years of experience, ACTEON is still the leader in piezo ultrasonics







- Compact and lightweight module for easy installation
- Choice of several ergonomic handpieces that are easily disconnectable and autoclavable
- Cords are light and flexible, with perfect cord/handpiece balance and control during the procedure
- Extremely precise electronic power settings for unbeatable power and performance
- Can be integrated into any delivery system or dental chair

ULTRASONICS



60x50x33mm / 75x49x33mm (LED) 28-36kHz 0.1 to 12W Supply Voltage: 24Vac / 30VA; 34Vdc Settings: analog or digital

SP NEWTRON LED MODULE



















NEWTRON STANDARD

Ø20 mm; 52g

ACCESS TO ALL ACTEON TIPS MISS









LED cold powerful light Ø 21mm; 73g











NEWTRON SLIM & B.LED

Non LED version available
Highly corrosion resistant (titanium made)
Ø 18mm; 48g
LED cold powerful light 65000 Lux





DIAGNOSTIC CAMERAS



SOPROCARE

Fluorescence technology highlights gingival inflammation, calculus and plaque

- PERIO mode: highlight plaque, calculus, and gingival inflammation
- CARIO mode: caries are detected as red, surrounding tissue is displayed in black and white
- **DAYLIGHT mode:** communicate more effectively with your patient and see details that are not visible with the naked eye







CARIO mode



DAYLIGHT mode




SOPROLIFE

Autofluorescence highlights decay & promotes minimally invasive treatment

- DIAGNOSTIC aid mode: identify the development of occlusal and interproximal carious lesions
- TREATMENT aid mode: perform minimally invasive treatment by preserving healthy tissue
- DAYLIGHT mode: from portrait to macrovision, obtain sharp images with the large depth of field



DIAGNOSTIC mode

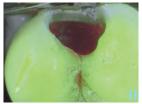

TREATMENT mode



DAYLIGHT mode



- Easy integration in all dental chairs
- Several configurations possible (USB, Video)
- Upgradable camera: one dock for all ACTEON cameras

INTRAORAL CAMERAS



SPEAK THE SAME LANGUAGE AS YOUR PATIENT



SOPR 0717

Macrovision reveals invisible: magnification up to 115 times

- Reveal micro fractures, infiltrations, lesions, everything that is not visible with the naked eye
- Extremely small camera head for easier access



Cracked tooth



Dental cavity preparation



Infiltrated occlusal



- Ease of use: point and shoot
- 105° angle of view for better exploration of distal areas
- Small, rounded shape for easy accessibility and maximum patient comfort
- Fixed focus with large depth of field, providing high-quality images



Intraoral



Intraoral



One tooth

INTRAORAL CAMERAS





- One digital USB 2.0 output
- 5VDC (from USB port)
- 2.5VA
- 48x48x30mm
- 22g













Several lengths available: 2m - 2.5m - 5m - 7m





STANDARD



Accessories



SHEATHS



SOPROTIPS

For SOPROLIFE and SOPROCARE




HANDPIECE HOLDER







A CAMERA FOR EVERY PRACTICE'S NEEDS

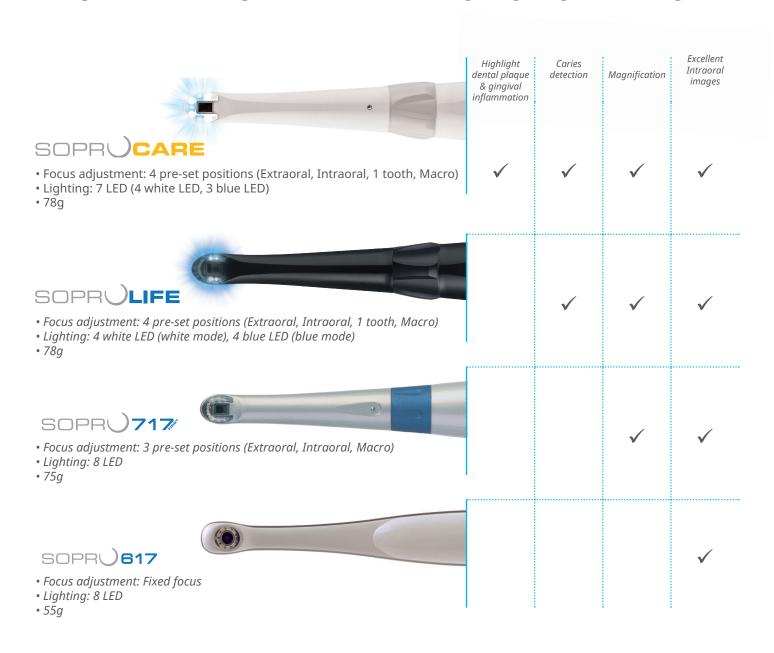

CURING LIGHTS




FAST CURING FOR LONG TERM ESTHETIC RESULTS




2 REMOVABLE AND AUTOCLAVABLE LIGHT GUIDES

- Ø 7.5mm: large restorations, posterior brackets, inlay cores
- Ø 5.5mm: lingual ortho, interproximal composites, small restorations



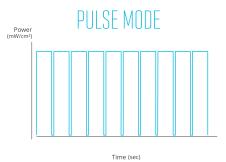

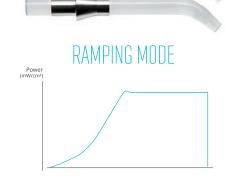


OUTPERFORMS IN POWER, EFFICIENCY & SPEED

- Robust aluminum handpiece
- Light and easy to grip
- Colored light shield for eye protection

High polymerization power

- Wide spectrum of wavelengths to cover the majority of dental photo-initiators
 - Standard and Supercharged 420-480nm
- Large power range (1250mW/cm² 3000mW/cm²)
- Focused energy and customizable modes for various polymerization depths

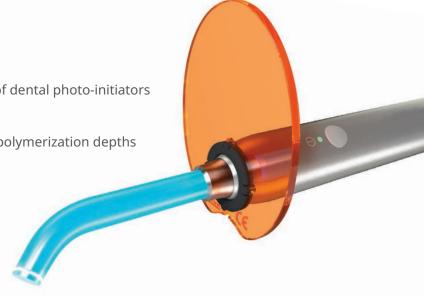



- All MiniLED™ use one single module with compact size, for easy integration and upgrades
- Various handpiece cords, shades and accessories



72x45x32mm Supply 24Vac / 20VA





Ø 4x6; 5x7mm Colors









MiniLED Standard Power: 1,250 - 2,000mW/cm² Spectrum: 420 - 480nm




MiniLED Supercharged Power: 2,000 - 3,000mW/cm² Spectrum: 420 - 480nm

AIR POLISHER



ALL-IN-ONE AIR POLISHER

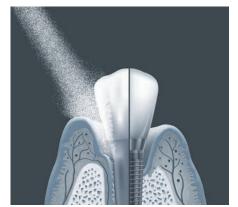

DUAL FUNCTIONS

- "SUPRA function" for supragingival polishing and stain removal
- "PERIO function" for sub-gingival polishing, periodontal maintenance and peri-implant treatments



4 NOZZLES FOR TREATMENT

AIR-N-GO Nozzles provide continuous, accurate and controlled spray. These nozzles offer power and efficiency for the removal of supra and subgingival biofilm.

- Sodium bicarbonate-based "CLASSIC" powders
- Calcium carbonate-based "PEARL" powder
- Glycine-based "PERIO" powder



AIRNGO

100% natural aromas or essential oils

A single handpiece to cover a range of care and treatment

- Fast start-up: direct connection to the high speed handpiece connector
- 360° swivel handpiece for fast and accurate treatment
- Parts can be dismantled for easy cleaning
- Quick-connect adapters to fit almost any handpiece adapter

HYGIENE



REDUCE THE RISK OF CROSS-CONTAMINATION

FIFXIBLE

Easy to bend for better access

ECONOMICAL

Saves time (no need for cleaning/disinfection/sterilization/traceability)

SAFE

- Reduced risk of direct or cross-contamination
- · No risk of ingestion: safety-locking device

DRY AIR

Two separate tubes guide air and water from the syringe to the end of the tip

UNIVERSAL

- Adapters** available for any type of equipment
- Adapters in anodized aluminum, propylux, stainless steel or resin, very easy to screw or clip onto any syringe










Sterilization of syringes is not completely effective!

Only disposable syringe tips like RISKONTROL® provide the best protection.

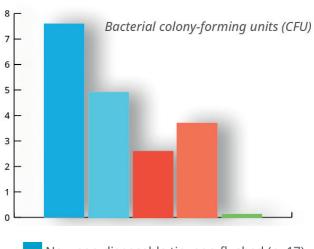

New non-disposable tip, non-flushed (n=17)

New non-disposable tip, flushed (n=17)

Used sterilizable tip, non-flushed (n=17)

Used sterilizable tip, flushed (n=17)

Riskontrol (n=17)

^{**} Class I medical devices

www.acteongroup.com

